


ITS.

Implants
trauma

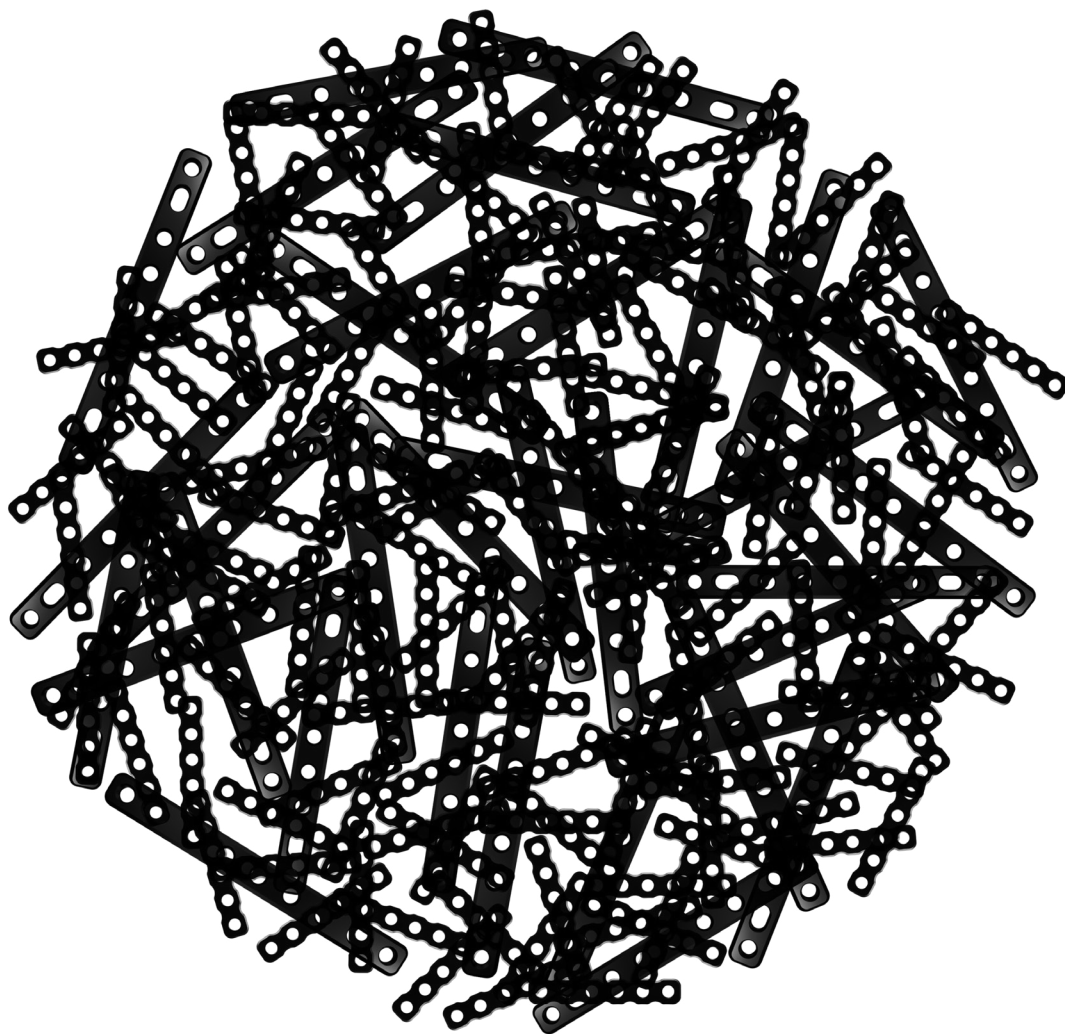

SLS

Straight Locking Plates System

CAUTION: Federal Law (USA) restricts this device to sale by or on the order of a board certified physician.

WARNING: If there is no sufficient bone healing, wrong or incomplete postoperative care, plate might break.

All ITS plates are preformed anatomically as a matter of principle. If adjustment of the plate to the shape of the bone is required, this is possible by carefully bending gently in one direction once. Particular care is required when bending in the region of a plate hole, as deformation of the plate may lead to a failure of the locking mechanism. The plate must not be buckled or bent several times. This is particularly important in the case of titanium implants, to prevent material fatigue and subsequent failure. The method of bending is the conscious responsibility of the operating doctor; I.T.S. GmbH can accept no liability whatsoever for this.

Contents

I. Introduction

- P. 5 Preface
- P. 6 Screws
- P. 7 Properties
- P. 8 Time of operation

2. Surgical Technique

- P. 8 Pre-operative patient preparation
- P. 8 Exposure
- P. 10 Straight Locking Plate 1.5mm
- P. 12 Straight Locking Plate 2.0mm
- P. 14 Straight Locking Plate 3.5/4.5mm
- P. 16 Postoperative treatment
- P. 16 Explantation

3. Information

- P. 17 Locking
- P. 17 Anodization Type II
- P. 18 Order list
- P. 22 Notes

Introduction




○ Preface

ITS. provides various Straight Locking Plates with different plate strengths and lengths covering a wide range of indications for shaft fractures.

The special feature of these implants is the free choice of screw placement. The user is able to set any desired screw in any hole either locking or non-locking screw (except long hole).

The free choice of screw angulation ($\pm 15^\circ$, see page 17) provides an advantage in fracture treatment, especially in the case of complex fractures.

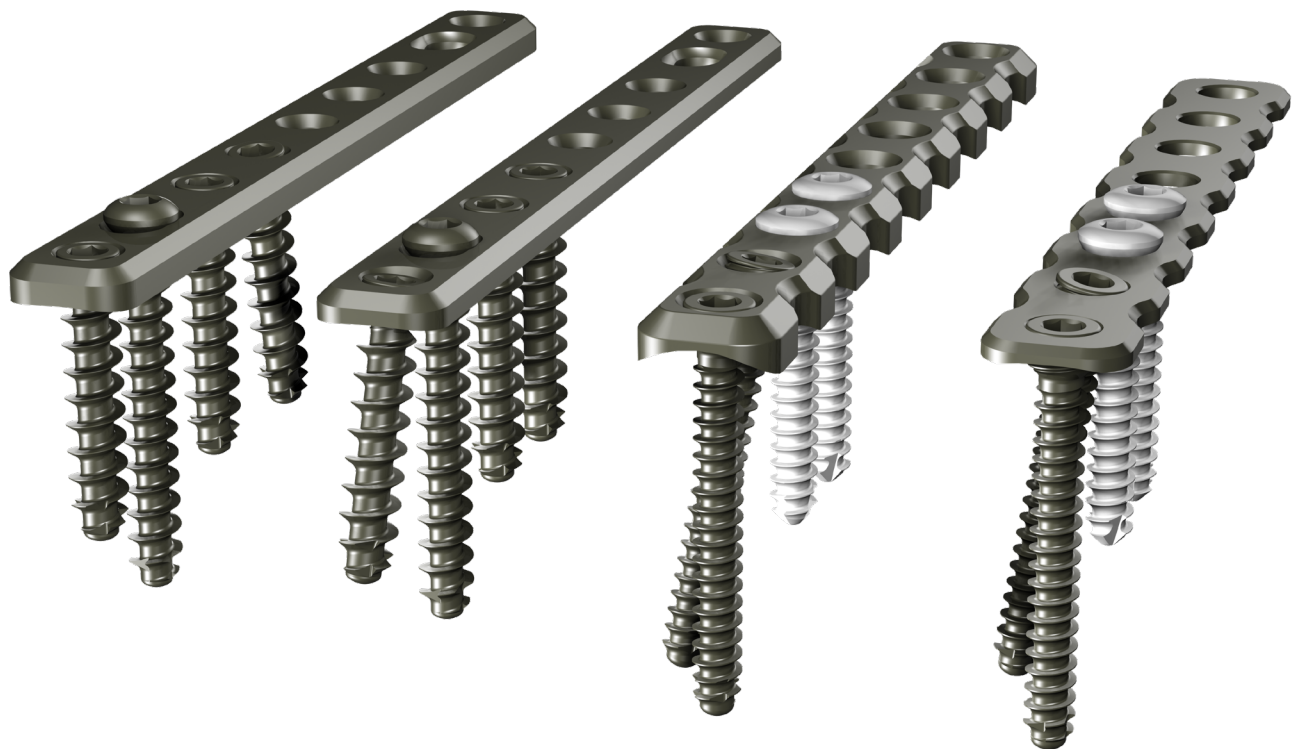


○ Screws Straight Locking Plate 1.5/2.0mm

3735I-XX-N Cortical Screw, locking, D=3.5mm, SH

61273-I00 Spiral Drill, D=2.7mm, L=100mm, AO Connector


56252 Screwdriver, WS 2.5, self-holding sleeve

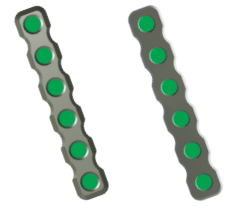

3235I-XX Cortical Screw, D=3.5mm

61273-I00 Spiral Drill, D=2.7mm, L=100mm, AO Connector



56252 Screwdriver, WS 2.5, self-holding sleeve

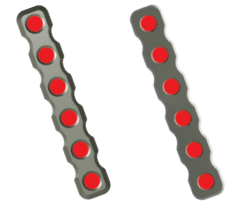

37422-XX-N Cancellous Screw, locking, D=4.2mm, SH

61253-II0 Spiral Drill, D=2.5mm, L=110mm, AO Connector

56252 Screwdriver, WS 2.5,
self-holding sleeve

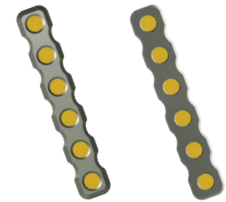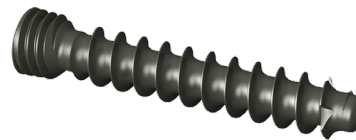


○ Screws Straight Locking Plate 3.5/4.5mm

37455-XX Cortical Screw, locking, D=4.5mm

61323-I45 Spiral Drill, D=3.2mm, L=145mm, AO Connector


56352-SH Screwdriver, WS 3.5, conic, self-holding

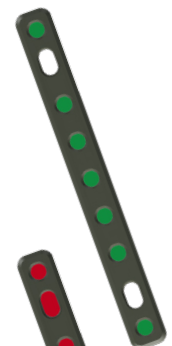

32455-XX Cortical Screw, D=4.5mm

61323-I45 Spiral Drill, D=3.2mm, L=145mm, AO Connector


56352-SH Screwdriver, WS 3.5, conic, self-holding

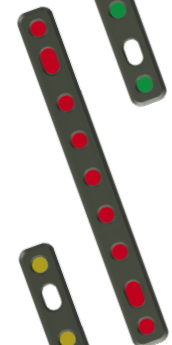

37592-XX Cancellous Screw, locking, D=5.9mm

61353-II0 Spiral Drill, D=3.5mm, L=110mm, AO Connector


56352-SH Screwdriver, WS 3.5, conic, self-holding

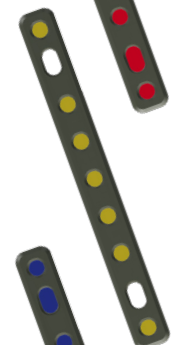

3059I-XX Cancellous Screw, D=5.9mm

61353-II0 Spiral Drill, D=3.5mm, L=110mm, AO Connector

56352-SH Screwdriver, WS 3.5, conic, self-holding

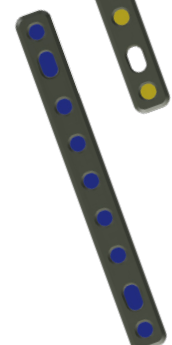


○ Properties

Properties of the material:

- Plate material: Titanium
- Material of screws: TiAl6V4 ELI
- Easier removal of the implant after the fracture has healed
- Improved fatigue strength of the implant
- Reduced risk of cold welding
- Reduced risk of inflammation and allergy

Properties of the implant:

- Multi-directional locking
- Anatomical plate design
- Accurate fitting of the plate on the bone
- Plate lengths, 1.5mm: 4, 5, 6, 7, 8, 9-hole
- Plate lengths, 2.0mm: 4, 5, 6, 7, 8, 9, 11, 13-hole
- Plate lengths, 3.5/4.5mm: 5, 6, 7, 8, 9, 10, 11, 12, 13-hole
- Straight Plate 3.5/4.5mm: Fracture gap compression up to 8mm

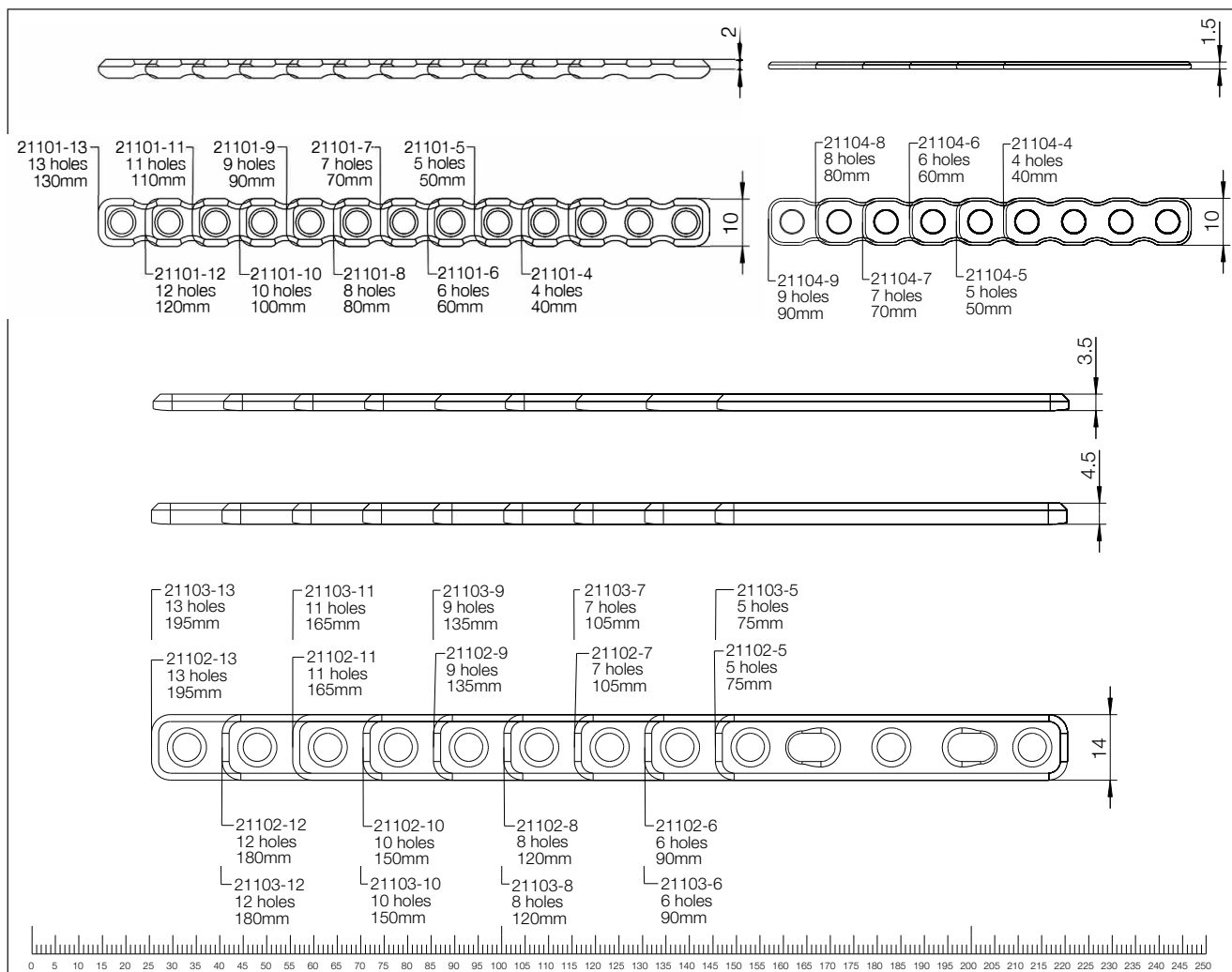

◦ Time of operation

- ♦ Immediately after trauma or delayed
- ♦ After regression of swelling

◦ Pre-operative patient preparation

- ♦ General anaesthesia, regional anaesthesia or combined can be used
- ♦ Tourniquet/partial deprivation of the blood supply

◦ Exposure

- ♦ Intra-operative x-ray fluoroscopy observation is required during the entire procedure
- ♦ Incision to the limb subchondral bone fracture site
- ♦ Transection of musculature if possible along the course of muscle fibers

2.

◦ Straight Locking Plate 1.5mm

Indications:

- ♦ Comminuted fractures
- ♦ Supracondylar fractures
- ♦ Intra-articular and extra-articular condylar fractures
- ♦ Fractures in osteopenic bone
- ♦ Non-unions, mal-unions
- ♦ Fracture or osteotomy of the tibia, fibula, femoral condyle, acetabulum, humerus, ulna, middle hand and middle foot bones
- ♦ Treatment of calcaneal
- ♦ Hip arthrodesis

Contraindications:

- ♦ All situations that do not allow osteosynthesis

◦ Plate insertion

- ♦ Anatomical reduction of the fracture segments
- ♦ Temporary fixation of the plate using guide wires
- ♦ Subsequent control under fluoroscopy
- ♦ Optionally, the plate can be stabilized using the ITS. Temporary Plate Holder (**58164-150**)



○ Placement of the screws

Use the drill guide, D=2.7/2.0mm (**62202**) to bore with the spiral drill, D=2.7mm, L=100mm, AO Connector (**61273-100**) into the plate holes closest to the fracture.

Use the screwdriver, WS 2.5, self-holding sleeve (**56252**) to insert D=3.5mm cortical screws (**3235I-XX**) of appropriate lengths determined previously with the depth gauge, solid small fragment screws (**59022**).

Attention: It is recommended that locking screws are not used close to the fracture.



Subsequently, place either D=3.5mm cortical or locking cortical screws (**3235I-XX/3735I-XX-N**) in plate holes far from the fracture.

Finally, control plate position under fluoroscopy.



◦ Straight Locking Plate 2.0mm

Indications:

- Comminuted fractures
- Supracondylar fractures
- Intra-articular and extra-articular condylar fractures
- Fractures in osteopenic bone
- Non-unions, mal-unions
- Fracture or osteotomy of the tibia, fibula, femoral condyle, acetabulum, humerus, ulna, middle hand and middle foot bones
- Treatment of calcaneal
- Hip arthrodesis

Contraindications:

- All situations that do not allow osteosynthesis

◦ Plate insertion

- Anatomical reduction of the fracture segments
- Temporary fixation of the plate using guide wires
- Subsequent control under fluoroscopy
- Optionally, the plate can be stabilized using the ITS. Temporary Plate Holder (**58164-150**)



○ Placement of the screws

Use the drill guide, D=2.7/2.0mm (**62202**) to bore with the spiral drill, D=2.7mm, L=100mm, AO Connector (**61273-100**) into the plate holes closest to the fracture.


Use the screwdriver, WS 2.5, self-holding sleeve (**56252**) to insert D=3.5mm cortical screws (**3235I-XX**) of appropriate lengths determined previously with the depth gauge, solid small fragment screws (**59022**).

Attention: It is recommended that locking screws are not used close to the fracture.



Subsequently, place either D=3.5mm cortical or locking cortical screws (**3235I-XX/3735I-XX-N**) in plate holes far from the fracture.

Finally, control plate position under fluoroscopy.

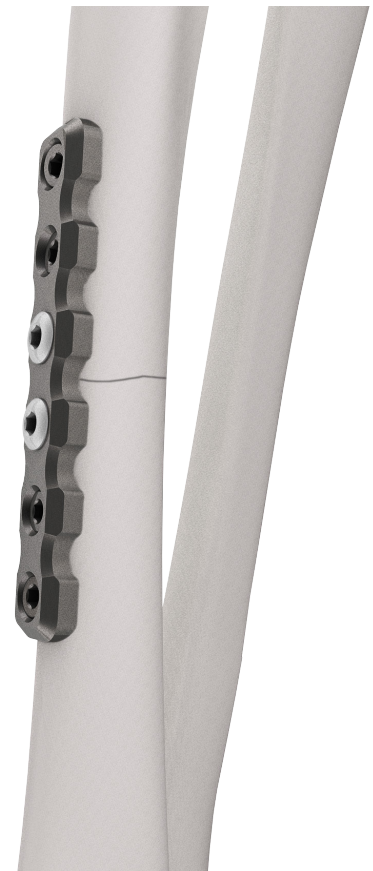

◦ Straight Locking Plate 3.5/4.5mm

Indications:


- Comminuted fractures
- Supracondylar fractures
- Intra-articular and extra-articular condylar fractures
- Fractures in osteopenic bone
- Non-unions, mal-unions
- Fracture or osteotomy of the tibia, fibula, femoral condyle, acetabulum, humerus, ulna, middle hand and middle foot bones
- Treatment of calcaneal
- Hip arthrodesis

Contraindications:

- All situations that do not allow osteosynthesis


◦ Plate insertion

- Anatomical reduction of the fracture parts
- Temporary fixation of the plate using forceps
- Subsequent control under fluoroscopy

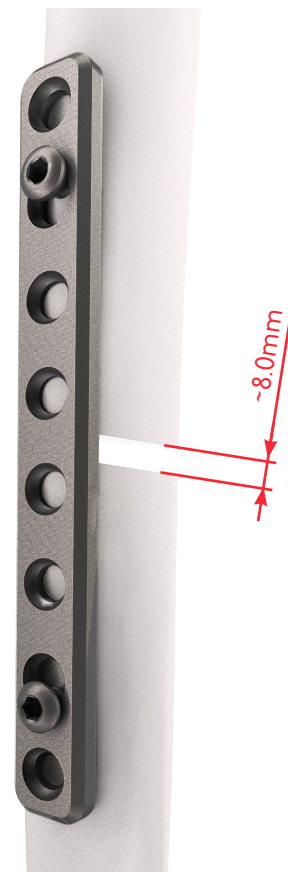

○ Placement of the screws


For compression up to 8mm the D=4.5mm cortical screws or D=5.9mm cancellous screws (32455-XX/30591-XX) are used at both sides of the fracture for compression.

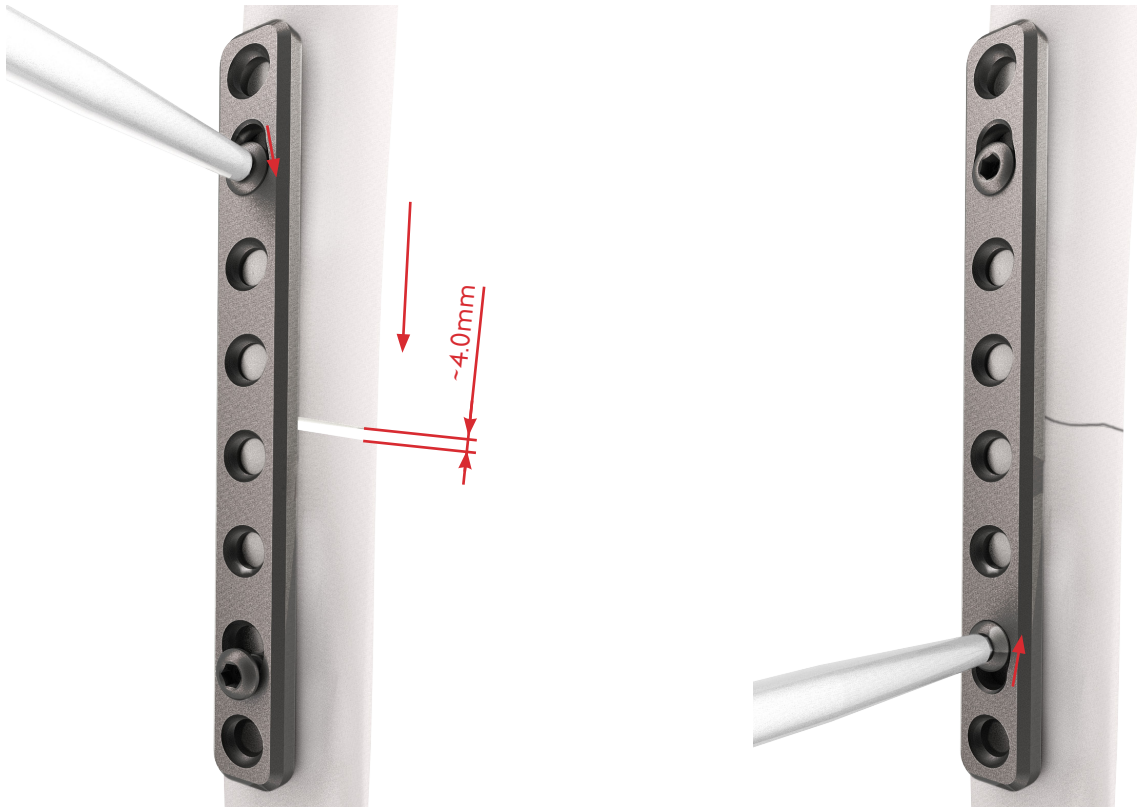

Subsequently, place screws in the remaining plate holes.

Attention: It is recommended that locking screws are not used close to the fracture.

Finally, control plate position under fluoroscopy.

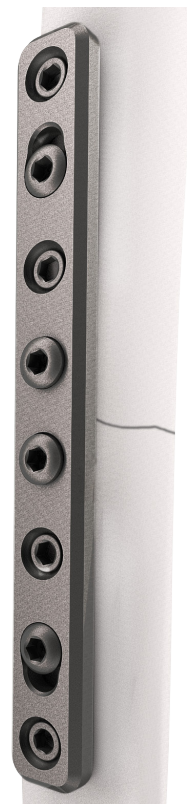

◦ Postoperative treatment

- Drainage is recommended for 12-24 hours to prevent postoperative hematoma
- Proper bandage dressing for 2 weeks (until the wound heals)
- Physical therapy for 5-7 weeks
- When a locking screw connection has been used, it is necessary to be aware that a diagnosis of non-union may be very delayed

◦ Explantation

If desired by the patient, the implant can be removed. Removal should be performed at the earliest 6 months - 1 1/2 years later or after radiographic verification of the healed bone.

The problem of cold welding was resolved by using a special surface treatment (for further information see page 17).


Information



○ Locking

Locking works because:

- Screw material (TiAlV) is slightly harder than plate material (Titanium Grade 2)
- Screw head **forms** thread into the plate (no cutting)

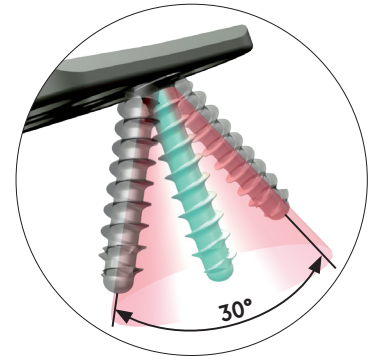

Benefits:

- $\pm 15^\circ$ and Locking
- No pre threading
- No cold welding
- No debris
- You can re-set the screw up to 3 times

○ Anodization Type II


Chemical process - anodization in a strong alkaline solution*

Type III anodization

- Layer thickness 60-200nm
- + Different colors
- Implant surface remains sensitive to:
Chipping
Peeling
Discoloration

Type II anodization

- Layer thickness 1000-2000nm
- + Film becomes an interstitial part of the titanium
- No visible cosmetic effect

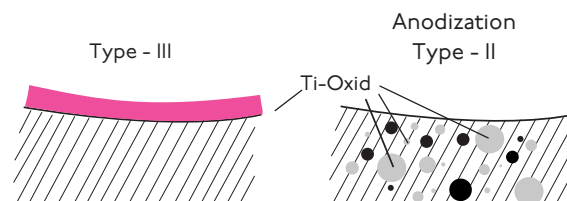

Anodization Type II leads to following benefits*

- Oxygen and silicon absorbing conversion layer
- Decrease in protein adsorption
- Closing of micro pores and micro cracks
- Reduced risk of inflammation and allergy
- Hardened titanium surface
- Reduced tendency of cold welding of titanium implants
- Increased fatigue resistance of implants
- Improved wear and friction characteristics

◦ Order list, Straight Locking Plate 1.5/2.0mm

Straight Plate, 1.5mm, 4-hole	21104-4
Straight Plate, 1.5mm, 5-hole	21104-5
Straight Plate, 1.5mm, 6-hole	21104-6
Straight Plate, 1.5mm, 7-hole	21104-7
Straight Plate, 1.5mm, 8-hole	21104-8
Straight Plate, 1.5mm, 9-hole	21104-9



Straight Plate, 2.0mm, 4-hole	21101-4
Straight Plate, 2.0mm, 5-hole	21101-5
Straight Plate, 2.0mm, 6-hole	21101-6
Straight Plate, 2.0mm, 7-hole	21101-7
Straight Plate, 2.0mm, 8-hole	21101-8
Straight Plate, 2.0mm, 9-hole	21101-9
Straight Plate, 2.0mm, 11-hole	21101-11
Straight Plate, 2.0mm, 13-hole	21101-13



Instruments and Screws, Straight Plate 1.5/2.0mm

Cortical Screw, D=3.5mm, L=10mm	32351-10
Cortical Screw, D=3.5mm, L=12mm	32351-12
Cortical Screw, D=3.5mm, L=14mm	32351-14
Cortical Screw, D=3.5mm, L=16mm	32351-16
Cortical Screw, D=3.5mm, L=18mm	32351-18
Cortical Screw, D=3.5mm, L=20mm	32351-20
Cortical Screw, D=3.5mm, L=22mm	32351-22
Cortical Screw, D=3.5mm, L=24mm	32351-24
Cortical Screw, D=3.5mm, L=26mm	32351-26
Cortical Screw, D=3.5mm, L=28mm	32351-28
Cortical Screw, D=3.5mm, L=30mm	32351-30
Cortical Screw, D=3.5mm, L=32mm	32351-32
Cortical Screw, D=3.5mm, L=34mm	32351-34
Cortical Screw, D=3.5mm, L=36mm	32351-36
Cortical Screw, D=3.5mm, L=38mm	32351-38
Cortical Screw, D=3.5mm, L=40mm	32351-40
Cortical Screw, D=3.5mm, L=42mm	32351-42
Cortical Screw, D=3.5mm, L=44mm	32351-44
Cortical Screw, D=3.5mm, L=46mm	32351-46
Cortical Screw, D=3.5mm, L=48mm	32351-48
Cortical Screw, D=3.5mm, L=50mm	32351-50
Cortical Screw, D=3.5mm, L=55mm	32351-55
Cortical Screw, D=3.5mm, L=60mm	32351-60



Cortical Screw, Locking, D=3.5mm, L=12mm, SH	37351-12-N
Cortical Screw, Locking, D=3.5mm, L=14mm, SH	37351-14-N
Cortical Screw, Locking, D=3.5mm, L=16mm, SH	37351-16-N
Cortical Screw, Locking, D=3.5mm, L=18mm, SH	37351-18-N
Cortical Screw, Locking, D=3.5mm, L=20mm, SH	37351-20-N
Cortical Screw, Locking, D=3.5mm, L=22mm, SH	37351-22-N
Cortical Screw, Locking, D=3.5mm, L=24mm, SH	37351-24-N
Cortical Screw, Locking, D=3.5mm, L=26mm, SH	37351-26-N
Cortical Screw, Locking, D=3.5mm, L=28mm, SH	37351-28-N



Cortical Screw, Locking, D=3.5mm, L=30mm, SH	37351-30-N
Cortical Screw, Locking, D=3.5mm, L=32mm, SH	37351-32-N
Cortical Screw, Locking, D=3.5mm, L=34mm, SH	37351-34-N
Cortical Screw, Locking, D=3.5mm, L=36mm, SH	37351-36-N
Cortical Screw, Locking, D=3.5mm, L=38mm, SH	37351-38-N
Cortical Screw, Locking, D=3.5mm, L=40mm, SH	37351-40-N

Cortical Screw, Locking, D=4.2mm, L=14mm, SH	37422-14-N
Cortical Screw, Locking, D=4.2mm, L=16mm, SH	37422-16-N
Cortical Screw, Locking, D=4.2mm, L=18mm, SH	37422-18-N
Cortical Screw, Locking, D=4.2mm, L=20mm, SH	37422-20-N
Cortical Screw, Locking, D=4.2mm, L=22mm, SH	37422-22-N
Cortical Screw, Locking, D=4.2mm, L=24mm, SH	37422-24-N
Cortical Screw, Locking, D=4.2mm, L=26mm, SH	37422-26-N
Cortical Screw, Locking, D=4.2mm, L=28mm, SH	37422-28-N
Cortical Screw, Locking, D=4.2mm, L=30mm, SH	37422-30-N
Cortical Screw, Locking, D=4.2mm, L=32mm, SH	37422-32-N
Cortical Screw, Locking, D=4.2mm, L=34mm, SH	37422-34-N
Cortical Screw, Locking, D=4.2mm, L=36mm, SH	37422-36-N
Cortical Screw, Locking, D=4.2mm, L=38mm, SH	37422-38-N
Cortical Screw, Locking, D=4.2mm, L=40mm, SH	37422-40-N
Cortical Screw, Locking, D=4.2mm, L=42mm, SH	37422-42-N
Cortical Screw, Locking, D=4.2mm, L=44mm, SH	37422-44-N
Cortical Screw, Locking, D=4.2mm, L=46mm, SH	37422-46-N
Cortical Screw, Locking, D=4.2mm, L=48mm, SH	37422-48-N
Cortical Screw, Locking, D=4.2mm, L=50mm, SH	37422-50-N
Cortical Screw, Locking, D=4.2mm, L=55mm, SH	37422-55-N
Cortical Screw, Locking, D=4.2mm, L=60mm, SH	37422-60-N



Screwdriver, WS 2.5, self-holding sleeve	56252
--	-------



Spiral Drill, D=2.7mm, L=100mm, AO Connector	61273-100
--	-----------



Depth Gauge, Solid Small Fragment Screws	59022
--	-------



Drill Guide, D=2.7/2.0mm	62202
--------------------------	-------



Temporary Plate Holder	58164-150
------------------------	-----------



Sterilization Tray, Straight Plate 1.5mm	50166
Sterilization Tray, Straight Plate 2.0mm	50166

For detailed cleaning and sterilization instructions, please refer to package insert.

◦ Order list, Straight Locking Plate 3.5/4.5mm

Straight Plate, 3.5mm, 5-hole	21102-5
Straight Plate, 3.5mm, 6-hole	21102-6
Straight Plate, 3.5mm, 7-hole	21102-7
Straight Plate, 3.5mm, 8-hole	21102-8
Straight Plate, 3.5mm, 9-hole	21102-9
Straight Plate, 3.5mm, 10-hole	21102-10
Straight Plate, 3.5mm, 11-hole	21102-11
Straight Plate, 3.5mm, 12-hole	21102-12
Straight Plate, 3.5mm, 13-hole	21102-13



Straight Plate, 4.5mm, 5-hole	21103-5
Straight Plate, 4.5mm, 6-hole	21103-6
Straight Plate, 4.5mm, 7-hole	21103-7
Straight Plate, 4.5mm, 8-hole	21103-8
Straight Plate, 4.5mm, 9-hole	21103-9
Straight Plate, 4.5mm, 10-hole	21103-10
Straight Plate, 4.5mm, 11-hole	21103-11
Straight Plate, 4.5mm, 12-hole	21103-12
Straight Plate, 4.5mm, 13-hole	21103-13



Cortical Screw, D=4.5mm, L=16mm	32455-16
Cortical Screw, D=4.5mm, L=20mm	32455-20
Cortical Screw, D=4.5mm, L=24mm	32455-24
Cortical Screw, D=4.5mm, L=28mm	32455-28
Cortical Screw, D=4.5mm, L=32mm	32455-32
Cortical Screw, D=4.5mm, L=36mm	32455-36
Cortical Screw, D=4.5mm, L=40mm	32455-40
Cortical Screw, D=4.5mm, L=44mm	32455-44
Cortical Screw, D=4.5mm, L=48mm	32455-48
Cortical Screw, D=4.5mm, L=52mm	32455-52
Cortical Screw, D=4.5mm, L=56mm	32455-56
Cortical Screw, D=4.5mm, L=60mm	32455-60
Cortical Screw, D=4.5mm, L=65mm	32455-65
Cortical Screw, D=4.5mm, L=70mm	32455-70
Cortical Screw, D=4.5mm, L=75mm	32455-75
Cortical Screw, D=4.5mm, L=80mm	32455-80
Cortical Screw, D=4.5mm, L=85mm	32455-85
Cortical Screw, D=4.5mm, L=90mm	32455-90



Cortical Screw, Locking, D=4.5mm, L=16mm	37455-16
Cortical Screw, Locking, D=4.5mm, L=20mm	37455-20
Cortical Screw, Locking, D=4.5mm, L=24mm	37455-24
Cortical Screw, Locking, D=4.5mm, L=28mm	37455-28
Cortical Screw, Locking, D=4.5mm, L=32mm	37455-32
Cortical Screw, Locking, D=4.5mm, L=36mm	37455-36
Cortical Screw, Locking, D=4.5mm, L=40mm	37455-40
Cortical Screw, Locking, D=4.5mm, L=44mm	37455-44
Cortical Screw, Locking, D=4.5mm, L=48mm	37455-48
Cortical Screw, Locking, D=4.5mm, L=52mm	37455-52
Cortical Screw, Locking, D=4.5mm, L=56mm	37455-56
Cortical Screw, Locking, D=4.5mm, L=60mm	37455-60
Cortical Screw, Locking, D=4.5mm, L=65mm	37455-65



Cortical Screw, Locking, D=4.5mm, L=70mm	37455-70
Cortical Screw, Locking, D=4.5mm, L=75mm	37455-75
Cortical Screw, Locking, D=4.5mm, L=80mm	37455-80
Cortical Screw, Locking, D=4.5mm, L=85mm	37455-85
Cortical Screw, Locking, D=4.5mm, L=90mm	37455-90

Cancellous Screw, Locking, D=5.9mm, L=16mm	37592-16
Cancellous Screw, Locking, D=5.9mm, L=20mm	37592-20
Cancellous Screw, Locking, D=5.9mm, L=24mm	37592-24
Cancellous Screw, Locking, D=5.9mm, L=28mm	37592-28
Cancellous Screw, Locking, D=5.9mm, L=32mm	37592-32
Cancellous Screw, Locking, D=5.9mm, L=36mm	37592-36
Cancellous Screw, Locking, D=5.9mm, L=40mm	37592-40
Cancellous Screw, Locking, D=5.9mm, L=44mm	37592-44
Cancellous Screw, Locking, D=5.9mm, L=48mm	37592-48
Cancellous Screw, Locking, D=5.9mm, L=52mm	37592-52
Cancellous Screw, Locking, D=5.9mm, L=56mm	37592-56
Cancellous Screw, Locking, D=5.9mm, L=60mm	37592-60
Cancellous Screw, Locking, D=5.9mm, L=65mm	37592-65
Cancellous Screw, Locking, D=5.9mm, L=70mm	37592-70
Cancellous Screw, Locking, D=5.9mm, L=75mm	37592-75
Cancellous Screw, Locking, D=5.9mm, L=80mm	37592-80
Cancellous Screw, Locking, D=5.9mm, L=85mm	37592-85
Cancellous Screw, Locking, D=5.9mm, L=90mm	37592-90



Cancellous Screw, D=5.9mm, L=16mm, Threaded	30591-16
Cancellous Screw, D=5.9mm, L=20mm, Threaded	30591-20
Cancellous Screw, D=5.9mm, L=24mm, Threaded	30591-24
Cancellous Screw, D=5.9mm, L=28mm, Threaded	30591-28
Cancellous Screw, D=5.9mm, L=32mm, Threaded	30591-32
Cancellous Screw, D=5.9mm, L=36mm, Threaded	30591-36
Cancellous Screw, D=5.9mm, L=40mm, Threaded	30591-40
Cancellous Screw, D=5.9mm, L=44mm, Threaded	30591-44
Cancellous Screw, D=5.9mm, L=48mm, Threaded	30591-48
Cancellous Screw, D=5.9mm, L=52mm, Threaded	30591-52
Cancellous Screw, D=5.9mm, L=56mm, Threaded	30591-56
Cancellous Screw, D=5.9mm, L=60mm, Threaded	30591-60
Cancellous Screw, D=5.9mm, L=65mm, Threaded	30591-65
Cancellous Screw, D=5.9mm, L=70mm, Threaded	30591-70
Cancellous Screw, D=5.9mm, L=75mm, Threaded	30591-75
Cancellous Screw, D=5.9mm, L=80mm, Threaded	30591-80
Cancellous Screw, D=5.9mm, L=85mm, Threaded	30591-85
Cancellous Screw, D=5.9mm, L=90mm, Threaded	30591-90





Screwdriver, WS 3.5, Conic, Self-holding	56352-SH
--	----------







Spiral Drill, D=3.2mm, L=145mm, AO Connector	61323-145
Spiral Drill, D=3.5mm, L=110mm, AO Connector	61353-110



For detailed cleaning and sterilization instructions, please refer to package insert.

Depth Gauge, Solid Small Fragment Screws	59022	
Drill Guide, D=2.5/3.5mm	62252	
Sterilization Tray, Straight Plate 3.5/4.5mm	50223	

Optional (on request)*

Cortical Screw, D=4.5mm, L=95mm	32455-95	
Cortical Screw, D=4.5mm, L=100mm	32455-100	
Cortical Screw, D=4.5mm, L=105mm	32455-105	
Cortical Screw, D=4.5mm, L=110mm	32455-110	
Cortical Screw, D=4.5mm, L=115mm	32455-115	
Cortical Screw, D=4.5mm, L=120mm	32455-120	
Cortical Screw, Locking, D=4.5mm, L=95mm	37455-95	
Cortical Screw, Locking, D=4.5mm, L=100mm	37455-100	
Cortical Screw, Locking, D=4.5mm, L=105mm	37455-105	
Cortical Screw, Locking, D=4.5mm, L=110mm	37455-110	
Cortical Screw, Locking, D=4.5mm, L=115mm	37455-115	
Cortical Screw, Locking, D=4.5mm, L=120mm	37455-120	
Cancellous Screw, Locking, D=5.9mm, L=95mm	37592-95	
Cancellous Screw, Locking, D=5.9mm, L=100mm	37592-100	
Cancellous Screw, Locking, D=5.9mm, L=105mm	37592-105	
Cancellous Screw, Locking, D=5.9mm, L=110mm	37592-110	
Cancellous Screw, Locking, D=5.9mm, L=115mm	37592-115	
Cancellous Screw, Locking, D=5.9mm, L=120mm	37592-120	
Cancellous Screw, D=5.9mm, L=95mm, Threaded	30591-95	
Cancellous Screw, D=5.9mm, L=100mm, Threaded	30591-100	
Cancellous Screw, D=5.9mm, L=105mm, Threaded	30591-105	
Cancellous Screw, D=5.9mm, L=110mm, Threaded	30591-110	
Cancellous Screw, D=5.9mm, L=115mm, Threaded	30591-115	
Cancellous Screw, D=5.9mm, L=120mm, Threaded	30591-120	

For detailed cleaning and sterilization instructions, please refer to package insert.

Notes



I.T.S. USA
1778 Park Avenue N, Suite 200
Maitland, FL 32751

Tel.: 877 - 971 - 8054
Fax: 877 - 971 - 8056
info@its-implantusa.com
www.its-implantusa.com



Order No. SLS-OP-0725-US
Edition: February/2025

© ITS. GmbH Graz/Austria 2025.
Subject to technical alterations, errors and misprints excepted.